

Permobil

F5 Corpus[®] VS



permobil

Stand without compromise

The F5 Corpus® VS power wheelchair offers Permobil power standing functionality on a front-wheel base for enhanced stability and ride comfort.

Stand confident

The F5 Corpus VS has a fully independent suspension on all wheels offering comfort, traction and stability, especially beneficial for active outdoor driving. Easy adjustment of the standing sequences

From face-to-face interactions to improved health outcomes, the F5 Corpus VS empowers you to stand with confidence, while our light, swing-away knee supports offer improved comfort and flexible on-chair storage for greater independence.

Stand in up to 80°

F5 Corpus VS has a standing angle of up to 80°, enabling you to independently perform activities of daily living, such as cooking, shopping and personal care.

Designed for independence

With a lightweight design, our swing-away knee supports can be moved and handled independently using the easy to reach release mechanism even with limited hand function. With two pad sizes, choose the option that best matches you and easily fine tune the pad placement as needed with the integrated BOA® dial* for optimal fit and comfort to your needs.

Programming simplified

QuickConfig is a new innovative app designed for service technicians, clinicians, and providers that offers a simpler way to configure parameters and an intuitive interface to set up F5 Corpus VS with your preferred driving and seating characteristics and standing sequence.



Base width
630 mm

Length incl.
Anti-tippers

*BOA® Fit System delivers micro-adjustable, precision fit in a simple to use dial

F5 Corpus[®] VS Specifications

Max. speed	6, 10 and 12 km/h
Range (G24: 85 Ah batteries)	40 km*
Max. user weight	136 kg
Obstacle climbing (with run up)	50 mm (75 mm)
Min. turning diameter	1490 mm
Turning in corridor	1200 mm
Suspension	Independent on all wheels
Max. standing angle	80°
R-net electronics Power Platform	120 A
Width incl. seat	650–790 mm
Height	1090–1170 mm
Min. transport length incl. Anti-tippers	965 mm
Min. transport height (folded backrest)	825 mm
Weight incl. batteries 85 Ah	206 kg
Weight of batteries 85 Ah	2x23 kg
Seat height (with electric seat elevator)	450–800 mm
Seat depth	370–570 mm (by 25 mm)
Seat width	420–570 mm (by 50 mm)
Backrest width	360/410/460/510 mm
Backrest height	470, 545–670 mm (by 25 mm)
Armrest height	215–330 mm
Distance between armrests	380–480/480–580 mm

*Driving conditions.

Measurements are approximating and subject to change without notice



Sit to Stand Seat Function



Lying Position to stand seat Function




Power Recline of 85-180°



Seat Elevator
Active Height 0-350 mm



Posterior Tilt 0-50°

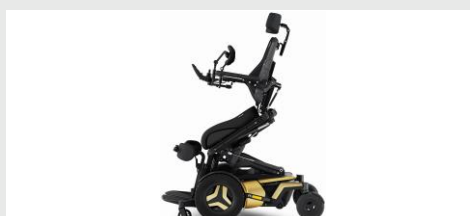

Anterior Tilt
Active Reach 5-45°

Personalise your ride



Standing supports
face-to face interactions
and enables improved health
outcomes.

Permobil

F5 Corpus[®] VS

Explore the vast array of standing wheelchair benefits with the F5 Corpus VS. It has all the performance features of the F5 Corpus combined with superior power standing functionality. The *Agile* suspension enables you to confidentially drive while standing at up to 5 km/h.

permobil



m.g.
MEDICOM SERVICES LTD

Tel.: 77771314, 24642070, Mob.: 99455233
e-mail: m.g.medicom@cytanet.com



GR3964-0.1UK230531