

Permobil

M3 Corpus[®]



permobil

More freedom

The M3 Corpus® has a compact base offering great maneuverability in tight spaces while with remarkable outdoor performance. It comes with comprehensive functionality and stability in various seating positions.

Corpus® – follows the contours of the body

The Corpus® seating system features an ergonomic design that follows the contours of the body. Improved comfort, pelvic positioning, lateral stability and immersion are all key objectives of the new Corpus cushions.

Corpus cushions feature softer dual-density foam on the backrest, a re-designed pelvic well and Stretch-Air™, a breathable, water-resistant cushion fabric that is machine-washable at 60°C.

Active Reach

Tilts the seat forward by up to 20° and raises it to extend functional reach, aid in transfers and your daily activities.

Active Height

Enables you to drive up to 5 km/h while fully elevated and combines seamlessly with *Active Reach* to give you the ultimate advantage in accessing your environment.

Power recline

Smooth and safe power recline with seamless sheer reduction for any user up to 150kg.

Programming simplified

QuickConfig is a new innovative app designed for service technicians, clinicians, and providers that offers a simpler way to configure parameters and an intuitive interface to set up M3 Corpus with your preferred driving and seating characteristics.



Active Reach truly putting the world within reach




The power recline with seamless sheer reduction moves between positions

M3 Corpus[®] Specifications

Max. speed	6/8/10 km/h
Range (G34: 65Ah / G24: 85Ah batteries)	30-40 km*
Max. user weight	150 kg
Obstacle climbing (with run up)	50 mm (75 mm)
Min. turning diameter	1120 mm
Turning in corridor	1100 mm
Suspension	<i>Smooth</i>
R-net electronics Power Platform	120 A
Width incl. seat	650-790 mm
Height	960-1170 mm
Min. transport length	855 mm
Min. transport height (folded backrest)	825 mm
Weight incl. batteries 65 Ah	176 kg
Weight of batteries 65 Ah	2x17 kg
Seat height (with electric seat elevator)	450-750 mm
Seat depth	370-570 mm (by 25 mm)
Seat width	420-570 mm (by 50 mm)
Backrest width	360/410/460/510 mm
Backrest height	470, 545-670 mm (by 25 mm)
Distance between armrests	380-480/480-580 mm

Personalise your ride

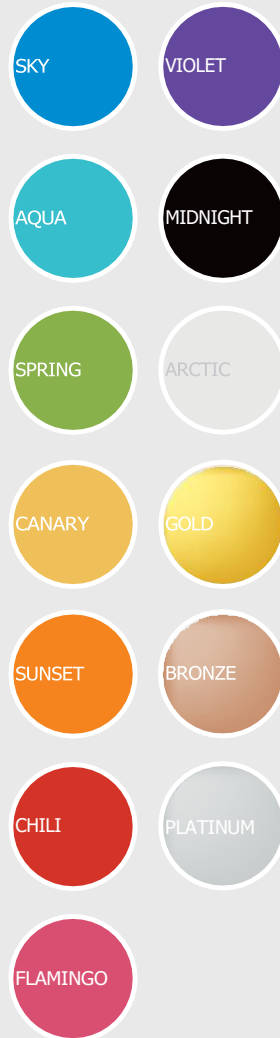

*Driving conditions and choice of battery type.
Measurements are approximates and subject to change without notice.

Explore the ultimate advantage to access the environment.

Permobil

M3 Corpus[®]

Meticulously engineered the M3 Corpus offers reliability, comfort, and great manoeuvrability in addition to the ultimate advantage to access the environment.

The M3 Corpus will simply offer more freedom.



permobil



m.g.
MEDICOM SERVICES LTD

Tel.: 77771314, 24642070, Mob.: 99455233
e-mail: m.g.medicom@cytanet.com



M3CORPUSDATASHEETUK230816