

INGLESE - cod.22.120 -20/09/2022



ORMESA®
• MADE IN ITALY SINCE 1980 •

DYNAMICO



ORMESA

DYNAMICO

Designed for your comfort

Dynamico by Ormesa is a gait trainer for adults and children, ideal for both outdoor and indoor use.

Dynamico can comfortably accommodate people with mild or moderate physical disability.

In the case of users with sufficient torso control, the chest support can be easily removed; for users with higher support needs **Ormesa has designed a wide range of additional components that make Dynamico highly adaptable.**

PRODUCT
MARKED 

DYNAMICO is a MEDICAL DEVICE that COMPLIES WITH THE SAFETY REQUIREMENTS defined with the REGULATION (EU) 2017/745



Built to last


A versatile aid with proven technology

Dynamico by Ormesa features a solid steel **frame** adjustable in height, **designed to withstand years of use and follow the user in all stages of the rehabilitation process.**

There are references on the steel tube for levelled height adjustment of pelvic and thoracic support and it has holes both for the application of weight bars (810) and for the adjustment of the abductor saddle (815).

Pelvic and thoracic supports, adjustable in height, help to make **Dynamico a valuable aid in multi-user environments** such as schools, clinics or nursing facilities.

Dynamico is available in **two different models:** Dynamico **for indoor use** and Dynamico **for outdoor use.**

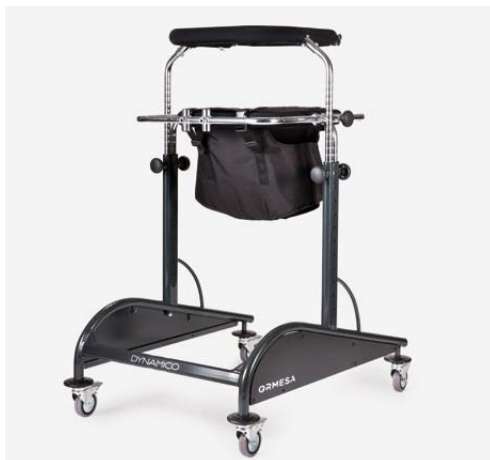

for indoor use available in 5 sizes

Indoor model is equipped with four pivoting wheels with independent brakes, which make the aid easy to manage inside the home, school or therapy spaces.



for outdoor use available in 3 sizes

Outdoor model can be used on roads and terrain without steep slopes in complete safety, thanks to its handy manual brakes.

Technical features and adjustments





Handwheels

For the adjustment in height of the thoracic and pelvic support.



Frame

Non-toxic epoxy varnished chromium-plated frame.



Hand brakes

Available only for Dynamico Outdoor use.



Harness

Supporting the user in case of fall. It's adjustable in height thanks to some belts and in girth with velcro strips. It can be easily removed for washing.



Bumpers

Semi-elastic PVC bumpers to protect from impacts.



Wheels

Indoor Dynamico

8cm, 4 castor wheels. They are all provided with independent brakes.

Outdoor Dynamico

17 cm, front pivoting wheels. 30 cm, back fixed wheels.



Handlebar

adjustable in height and forward-backward.



Thoracic support

padded support for the trunk, adjustable in height. It can be removed for people with sufficient trunk control. It allows an easy access thanks to back opening.



Pelvic support

padded support for the pelvis. It's adjustable either in height or in forward-backward direction in order to define the walking, moving of the barycentre forwards. It allows an easy access thanks to back opening.



Components



809 Arm supports

With depth, inclination, rotation, height and prone supination adjustment. Alternative to 813 Handlebar and 824 Tray.




943 Arm straps

With circumference adjustment.



824 Transparent tray

With circular recess. With height and forward-backward adjustment. Alternative to 809 Arm support and 813 Handlebar. Only for Dynamico Indoor use.

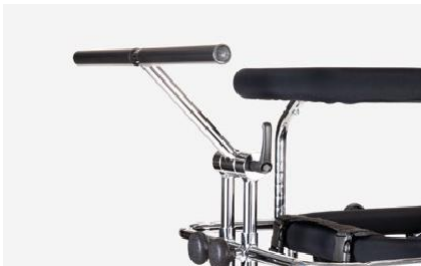

813 Adjustable handlebar

With height and forward-backward adjustment. Alternative to 809 Arm support and 824 Tray.



815 Abduction saddle

Adjustable in height through holes and in forward-backward. Available for sizes 2,3 and 4. Only for Indoor use.



814 Leg divider

Only for Dynamico Indoor use.

814E Leg divider

Only for Dynamico Outdoor use.

Components



811 Adjustable brake for front wheels

Helps the adjustment of flowability of the walking frame. Only for Indoor use.



812 Directional locks for back wheels

Helps the direction of the walking frame. Only for Indoor use.




810 Weight bars

Available for sizes 2, 3, 4 and 5.
Size 2 + 2,5 kg. x 2
Size 3 + 3,1 kg. x 2
Size 4 + 3,7 kg. x 2
Size 5 + 4,35 kg. x 2



Indoor Dynamico sizes and weight


	size 1	size 2	size 3	size 4	size 5
A: height ground to pelvic support	from 44 to 60 cm <i>from 17,3 to 23,6 in</i>	from 66 to 82 cm <i>from 25,9 to 32,2 in</i>	from 77 to 102 cm <i>from 30,3 to 40,1 in</i>	from 85 to 112 cm <i>from 33,4 to 44 in</i>	from 97 to 125 cm <i>from 38,1 to 49,2 in</i>
B: height ground to thoracic support	from 52 to 70 cm <i>from 20,4 to 27,5 in</i>	from 77 to 95 cm <i>from 30,3 to 37,4 in</i>	from 89 to 116 cm <i>from 35 to 45,6 in</i>	from 92 to 126 cm <i>from 36,2 to 49,6 in</i>	from 105 to 138 cm <i>from 41,3 to 54,3 in</i>
C: depth of inner supports	18 cm ~ 7 in	21 cm ~ 8,2 in	25 cm ~ 9,8 in	30 cm ~ 11,8 in	32 cm ~ 12,5 in TS* 41 cm ~ 16,1 in PS*
D: width of inner supports	25 cm ~ 9,8 in	28 cm ~ 11 in	30 cm ~ 11,8 in	35 cm ~ 13,7 in	41 cm ~ 16,1 in
E: circumference of supports	65 cm ~ 25,5 in	75 cm ~ 29,5 in	85 cm ~ 33,4 in	110 cm ~ 43,3 in	120 cm ~ 47,2 in
F: length encumbrance	79 cm ~ 31,1 in	84 cm ~ 33 in	88 cm ~ 34,6 in	95 cm ~ 37,4 in	106 cm ~ 41,7 in
G: overall width	66 cm ~ 25,9 in	66 cm ~ 25,9 in	66 cm ~ 25,9 in	71 cm ~ 27,9 in	77 cm ~ 30,3 in
H: user height	from 75 to 100 cm <i>from 29,5 to 39,3 in</i>	from 95 to 130 cm <i>from 37,4 to 51,1 in</i>	from 120 to 155 cm <i>from 47,2 to 61 in</i>	from 140 to 175 cm <i>from 55,1 to 68,8 in</i>	from 160 to 195 cm <i>from 62,9 to 76,7 in</i>
Total weight	12 kg ~ 26,4 lb	14,6 kg ~ 32,1 lb	15,6 kg ~ 34,3 lb	17 kg ~ 37,4 lb	19,7 kg ~ 43,4 lb
Maximum load	55 kg ~ 121,2 lb	70 kg ~ 154,3 lb	85 kg ~ 187,3 lb	100 kg ~ 220,4 lb	115 kg ~ 253,5 lb
Turning width	97 cm ~ 38,1 in	101 cm ~ 39,7 in	104 cm ~ 40,9 in	117 cm ~ 46 in	127 cm ~ 50 in

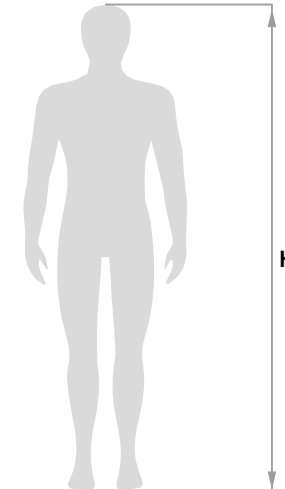

*TS: thoracic supp.
*PS: pelvic supp.

OutdoorDynamico sizes and weight

	size 3	size 4	size 5
A: height ground to pelvic support	from 85 to 110 cm <i>from 33,4 to 43,3 in</i>	from 93 to 120 cm <i>from 36,6 to 47,2 in</i>	from 106 to 133 cm <i>from 41,7 to 52,3 in</i>
B: height ground to thoracic support	from 97 to 124 cm <i>from 38,1 to 48,8 in</i>	from 100 to 134 cm <i>from 39,3 to 52,7 in</i>	from 116 to 146 cm <i>from 45,6 to 57,4 in</i>
C: depth of inner supports	25 cm ~ 9,8 in	30 cm ~ 11,8 in	32 cm ~ 12,5 in TS* 41 cm ~ 16,1 in PS*
D: width of inner supports	30 cm ~ 11,8 in	35 cm ~ 13,7 in	41 cm ~ 16,1 in
E: circumference of supports	85 cm ~ 33,4 in	110 cm ~ 43,3 in	120 cm ~ 47,2 in
F: length encumbrance	94 cm ~ 37 in	101 cm ~ 39,7 in	111 cm ~ 43,7 in
G: overall width	74 cm ~ 29,1 in	78 cm ~ 30,7 in	82 cm ~ 32,2 in
H: user height	from 120 to 155 cm <i>from 47,2 to 61 in</i>	from 140 to 175 cm <i>from 55,1 to 68,8 in</i>	from 160 to 195 cm <i>from 62,9 to 76,7 in</i>
Total weight	19,3 kg ~ 42,5 lb	20,4 kg ~ 44,9 lb	22,7 kg ~ 50 lb
Maximum load	85 kg ~ 187,3 lb	100 kg ~ 220,4 lb	115 kg ~ 253,5 lb

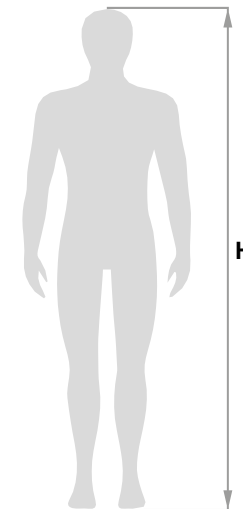

*TS: thoracic supp.
*PS: pelvic supp.

DYNAMICO





ORMESA

ORMESA[®]

• MADE IN ITALY SINCE 1980 •




m.g.

MEDICOM SERVICES LTD

Τηλ.: 77771314, 24642070, Κινητό: 99455233

www.medicomstore.com

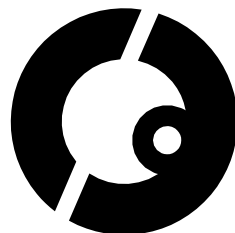

 complies with MDR (EU) 2017/745

COMPANY WITH
QUALITY SYSTEM
CERTIFIED BY DNV GL
= ISO 13485 =

Ormesa srl - Via delle Industrie 6/8 | 06034 Foligno (PG) - ITALY
Tel. +39 0742 22927 | FAX +39 0742 22637 | info@ormesa.com

www.ormesa.com